

# Thermo Scientific Synchronis Silica HPLC Columns

Consistent, predictable separations,  
Column after column, time after time

- Highly pure, high surface area silica
- Excellent reproducibility for normal phase chromatography
- Rigorously tested to ensure quality



## Specifications


<b>Particle size</b>	1.7 $\mu\text{m}$ , 3 $\mu\text{m}$ , 5 $\mu\text{m}$	<b>Carbon load</b>	–
<b>Pore size</b>	100 Å	<b>Endcapped</b>	–
<b>Surface area</b>	320 $\text{m}^2/\text{g}$	<b>USP classification</b>	L3
<b>pH range</b>	2-8		

## High efficiency columns for normal phase chromatography

Unbonded, highly pure base deactivated silica media that is the backbone of the Thermo Scientific™ Synchronis™ range of columns. With a high surface area, Synchronis Silica HPLC columns

are a powerful and efficient tool for the chromatography of non-polar and moderately polar organic compounds by normal phase chromatography. For high speed, high efficiency separations, Synchronis Silica HPLC columns are available with 1.7  $\mu\text{m}$  particle size.

## Application: Cetirizine (USP)

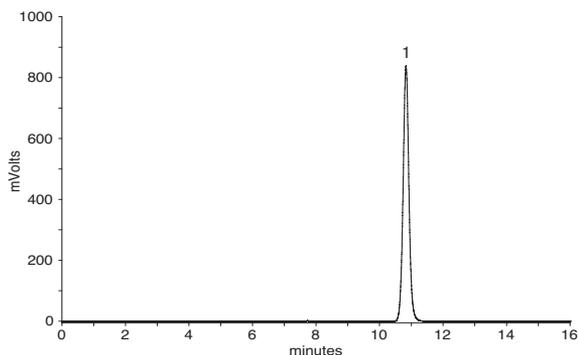


**Column:** Synchronis Silica 5  $\mu\text{m}$ , 250 mm x 4.6 mm  
**Mobile phase:** acetonitrile:water:sulfuric acid (93:6.6:0.4)  
**Flow rate:** 1.0 mL/min  
**Temperature:** 30 °C  
**Detection:** 230 nm  
**Injection volume:** 10  $\mu\text{L}$

1. Cetirizine

Parameter	USP Specification	Measured (6 replicate injections)
Tailing factor	< 2.0	1.05
%RSD peak area	< 2%	0.17%

**Thermo**  
SCIENTIFIC

## Normal phase retention and selectivity

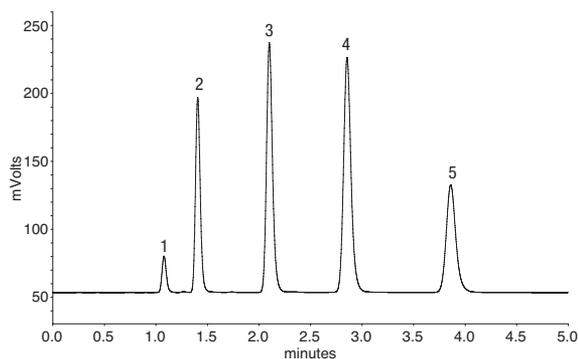

**Column:** Synchronis Silica 5  $\mu\text{m}$ , 150 mm x 4.0 mm  
**Mobile phase:** iso-octane:ethanol:water (85:14.7:0.3)  
**Flow rate:** 1.25 mL/min  
**Temperature:** 40 °C  
**Detection:** 254 nm  
**Injection volume:** 2  $\mu\text{L}$

1. Toluene  
 2. Nitrobenzene  
 3. o-nitroaniline  
 4. m-nitroaniline  
 5. p-nitroaniline

## Ordering Information

Description	Particle Size	Length (mm)	2.1 mm ID	3.0 mm ID	4.0 mm ID	4.6 mm ID
Synchronis Silica	1.7 $\mu\text{m}$	30	97002-032130	–	–	–
		50	97002-052130	97002-053030	–	97002-054630
		100	97002-102130	97002-103030	–	–
	3 $\mu\text{m}$	50	97003-052130	–	–	–
		100	–	–	–	97003-104630
		150	–	–	–	97003-154630
	5 $\mu\text{m}$	30	97005-032130	97005-033030	97005-034030	97005-034630
		50	97005-052130	97005-053030	97005-054030	97005-054630
		100	97005-102130	97005-103030	97005-104030	97005-104630
		150	97005-152130	97005-153030	97005-154030	97005-154630
250		97005-252130	97005-253030	97005-254030	97005-254630	
Drop-in guard cartridges (4/pk)	5 $\mu\text{m}$	10	97005-012101	97005-013001	97005-014001	–
UNIGUARD Guard Cartridge Holder			852-00	852-00	850-00	–

## Consistent, predictable separations: Column after column, time after time

Synchronis HPLC columns are manufactured, packed and tested in ISO 9001:2008 accredited facilities. Each lot of silica is tested for the physical properties of the silica support and only released for production if it meets the stringent test specifications.

Each bonded lot of chromatographic packing material is rigorously tested for primary and secondary interactions with the bonded phase.

New, enhanced, automated packing methods drive consistency even further and every column is individually tested to ensure that it meets the required quality.

These extensive testing and quality control procedures ensure the delivery of a consistent product, column after column.

[thermoscientific.com/Synchronis](http://thermoscientific.com/Synchronis)

© 2012 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

**USA and Canada** +1 800 332 3331  
**France** +33 (0)1 60 92 48 34  
**Germany** +49 (0) 2423 9431 20 or 21  
**United Kingdom** +44 (0)1928 534110  
**Japan** +81 3 5826 1615

**China** +86 21 68654588 +86 10 84193588  
 +86 20 83145199 800 810 5118  
**India** +91 22 6742 9494  
**Australia** 1 300 735 292 (free call domestic)  
**New Zealand** 0800 933 966 (free call domestic)  
**All Other Enquiries** +44 (0) 1928 534 050

**Technical Support**  
**North America** +1 800 332 3331  
**Outside North America** +44 (0) 1928 534 440

**Thermo**  
 SCIENTIFIC

Part of Thermo Fisher Scientific